



Pulmodyne BiTrac®
Shields & Masks



Proudly Providing Airway
Solutions for Over 30 years

Fixed and Select Options to Suit a Variety of Patient Needs

- The BiTrac Full Face and Nasal Masks come in a variety of sizes. Both product lines are available with fixed elbows in place
- Our BiTrac Select™ Full Face Masks come with interchangeable elbows and a variety of cushion styles, to suit multiple patient needs

BiTrac™ Nasal with fixed Anti-Asphyxia Elbow

Order No.	Description	Qty
313-9013W	Adult Small	10/CS
313-9014W	Adult Medium	10/CS
313-9015W	Adult Large	10/CS
313-9016W	Adult Extra Large	10/CS

BiTrac™ Full Face with fixed Anti-Asphyxia Elbow

313-9001W	Adult Small	10/CS
313-9002W	Adult Medium	10/CS
313-9003W	Adult Large	10/CS
313-9004W	Adult Extra Large	10/CS

BiTrac™ Full Face with Fixed Standard Elbow

313-9009W	Adult Small	10/CS
313-9010W	Adult Medium	10/CS
313-9011W	Adult Large	10/CS
313-9012W	Adult Extra Large	10/CS

BiTrac Select™ Full Face with Interchangeable Anti-Asphyxia with leak Elbow

Order No.	Description	Qty
313-9511	Adult Small	10/CS
313-9512	Adult Medium	10/CS
313-9513	Adult Large	10/CS
313-9514	Adult Extra Large	10/CS

BiTrac Select™ Full Face with Interchangeable Standard Elbow

313-9506	Adult Small	10/CS
313-9507	Adult Medium	10/CS
313-9508	Adult Large	10/CS
313-9509	Adult Extra Large	10/CS

BiTrac Select™ Full Face with Interchangeable Anti-Asphyxia Elbow

313-9501	Adult Small	10/CS
313-9502	Adult Medium	10/CS
313-9503	Adult Large	10/CS
313-9504	Adult Extra Large	10/CS
313-9593	ICE Anti-Asphyxia w/ Leak Elbow	10/CS
313-9594	ICE Anti-Asphyxia Nebulizer Elbow	10/CS
313-9595	ICE Standard Nebulizer Elbow	10/CS
313-9596	ICE Anti-Asphyxia Bronchoscopy Elbow	10/CS
313-9597	ICE Standard Bronchoscopy Elbow	10/CS
313-9598	ICE Anti-Asphyxia Elbow	10/CS
313-9599	ICE Standard Elbow	10/CS



Fixed and Select Options to Suit a Variety of Patient Needs

The BiTrac® MaxShield offers a total facial perimeter cushion made with a soft silicone. The total face application helps prevent skin breakdown on the bridge of the nose. Interchanging the BiTrac® MaxShield with the BiTrac® Full Face Mask will support pressure ulcer prevention.

- MaxShield comes in a wide variety of sizes and elbow options
- Assists in management of claustrophobia
- MaxShield comes in a Fixed line and a Select line with interchangeable elbows
- A quick and easy fit for patients



BiTrac MaxShield Select™ with Interchangeable Anti-Asphyxia Elbow

Order No.	Description	Qty
313-9559	Adult Extra Large	10/CS
313-9558	Adult Large	10/CS
313-9557	Adult Small	10/CS

BiTrac MaxShield Select™ with Interchangeable Anti-Asphyxia with Leak Elbow

313-9555	Adult Extra Large	10/CS
313-9552	Adult Large	10/CS
313-9554	Adult Small	10/CS

BiTrac MaxShield SE Select™ with Interchangeable Standard Elbow

313-9556	Adult Extra Large	10/CS
313-9553	Adult Large	10/CS
313-9551	Adult Small	10/CS

Interchangeable Elbows

Order No.	Description	Qty
313-9599	Standard Elbow (dual limb circuit)	10/CS
313-9593	Anti-Asphyxia with Leak Elbow (single limb circuit)	10/CS
313-9597	Standard Bronchoscopy Elbow (dual limb circuit)	10/CS

BiTrac MaxShield Select™ with fixed Anti-Asphyxia Elbow

313-9057W	Adult Extra Large	10/CS
313-9052W	Adult Large	10/CS
313-9054W	Adult Small	10/CS

BiTrac SE MaxShield™ with fixed Standard Elbow


313-9056W	Adult Extra Large	10/CS
313-9053W	Adult Large	10/CS
313-9055W	Adult Small	10/CS



Expanded MaxShield Sizing for Pediatric Patients

The MaxShield now comes in more sizes, XXS, XS, and Adult Small, to accommodate the pediatric patient population. Like the adult line, the pediatric line offers a full perimeter fit, which aids in reducing claustrophobia.

The elbows are interchangeable, with a child proof bayonet locking system to hold them in place. The Standard Elbow can be used with a dual limb circuit, and the Anti-Asphyxia with Leak Elbow can be used with a single limb circuit. The MaxShield also comes with a bonnet style headgear with a soft foam for patient comfort.

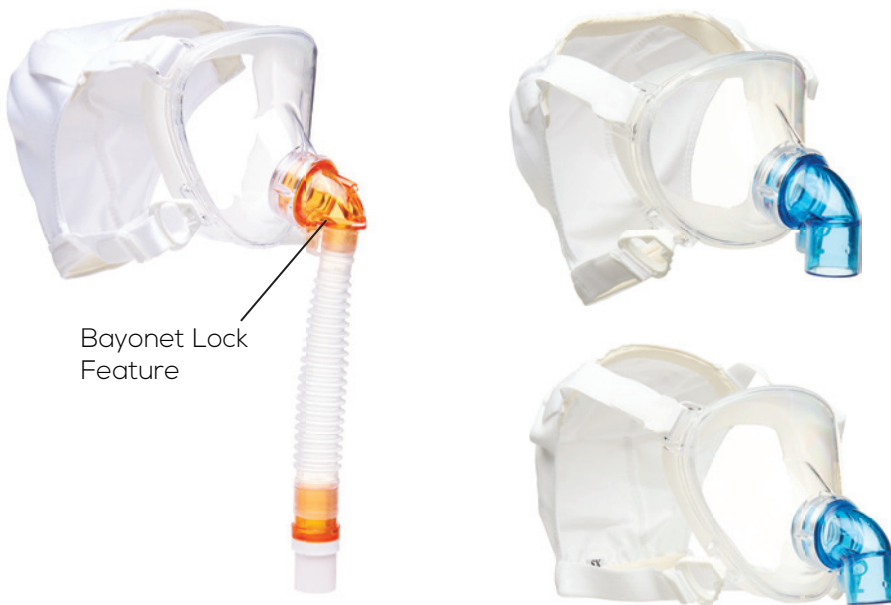

BiTrac® MaxShield Select™ Pediatric NIV Mask

Order No.	Size	Headgear	Elbow	Qty
313-9524	XXS	CAP Bonnet	Pediatric Standard	10/CS
313-9525	XS	CAP Bonnet	Pediatric Standard	10/CS
313-9526	XXS	CAP Bonnet	Pediatric Anti-Asphyxia with Leak	10/CS
313-9527	XS	CAP Bonnet	Pediatric Anti-Asphyxia with Leak	10/CS
313-9528	XXS	CAP Bonnet (Pediatric MaxShield Headgear Only)	-	10/CS
313-9529	XS	CAP Bonnet (Pediatric MaxShield Headgear Only)	-	10/CS
313-9602	-	None	Pediatric Anti-Asphyxia with Leak	10/CS



Set It and Concentrate on Caring for Your Patient

The GO-PAP is a disposable emergency CPAP device that runs off the barbed valve outlet on an O2 source. It operates at 10LPM and delivers approximately 90LPM at 30% FiO2. There is an adjustable PEEP(5/7.5/10cm H2O) valve integrated into the generator. The PEEP is independent of the flow, so you never have to worry about adjusting your flow to maintain your PEEP levels. Additionally, this completely disposable system offers integrated nebulization while maintaining a consistent flow to the patient.

- Runs at 10 LPM of oxygen, and generates up to 90 LPM of flow
- Has an approximate 30% FiO2
- Includes 3 options for PEEP of 5, 7.5, and 10 cm H2O
- Completely disposable system
- Run time of approximately 40 minutes
- Integrated nebulization without interrupting the delivered flow to the patient



GO-PAP™ System with BiTrac™ ED Mask and Neb-Connect™- Nebulizer Canister Included

Order No.	Headgear	Size	Qty
313-4604NC	CAP Adult	Medium	10/CS
313-4605NC	CAP Adult	Large	10/CS
313-4602NC	Standard Adult	Medium	10/CS
313-4603NC	Standard Adult	Large	10/CS

Neb-Connect™

Order No.	Recommended Flow Rate	Qty.
313-011111	6-8 LPM	30/CS

Neb-Connect™

Allows You to CPAP & Nebulize in ONE Tank!



GO-PAP™ System with BiTrac™ ED Mask and Neb-Connect™ - Nebulizer Canister Included

Order No.	Headgear	Size	Qty
313-4601N	Standard Adult	Small	10/CS
313-4602N	Standard Adult	Medium	10/CS
313-4603N	Standard Adult	Large	10/CS

GO-PAP™ System with BiTrac™ ED Mask and Neb-Connect™ - Nebulizer Canister Not Included

313-4601	Standard Adult	Small	10/CS
313-4602	Standard Adult	Medium	10/CS
313-4603	Standard Adult	Large	10/CS

